

Instruction Manual

e-pill[®] MedTime Station Automatic Pill Dispenser with Station and Cup

SKU 997025 | US Patent 9,198,834 B2 | US Patent 9,539,177 B2



Carefully review and save this setup guide



Medication Reminders

Station

(metal base)

To MOUNT the Station to a surface (Required)

Peel the plastic covering off the four sticky feet located on the bottom of the Station, then place the Station in your chosen area

To REMOVE the round pill dispenser from the metal Station

TURN the pill dispenser counter-clockwise until the e-pill logo is at 10 o'clock, then **LIFT** up

To SECURE the round pill dispenser to the metal Station

ALIGN the pill dispenser with the e-pill logo at 10 o'clock
TURN the pill dispenser clockwise until the e-pill logo is at 12 o'clock

e-pill® MedTime Station Locked Automatic Pill Dispenser with Station



Introduction

The e-pill MedTime Station is designed to remind and dispense pills according to the patient's dosing schedule; reminding by means of auditory alarm and flashing light. The dispenser will make the correct dose available at the correct time, keeping all other pills secure via a locked lid. The metal Station makes dispensing pills for those with dexterity concerns easy. No more spilled pills.

e-pill MedTime Station is locked but it is not tamperproof, therefore it should NOT be used to administer pain medications or Schedule II/III drugs.

Do not place e-pill MedTime Station in direct sunlight as it may cause operational malfunction.

How it Works

At dose time the alarm will sound, and the internal tray will rotate making the correct dose available. Pull the blue handle towards you to release the correct dose and silence the alarm.


Components



Locked Automatic Pill Dispenser (1)
Item: Not available separately



Medication Tray (1)
Item: 996030

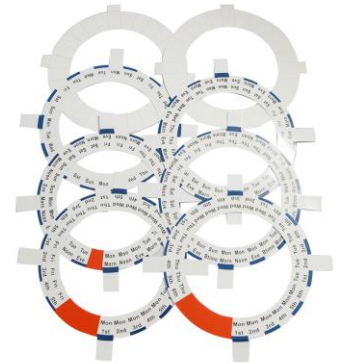

Dosage Template Set (1)
Item: 996035



Station (1)
Item: 991024



Stainless Steel Dose Cup (1)
Item: 992028



Keys (2)
Item: 993026



AA Batteries (4)
Item: Not available separately

Note: We recommend that you always keep e-pill MedTime Station locked. Two keys are included: One key for the caregiver and one backup key.

SETUP - STEP 1



Display and Buttons:

1. ON: Button 3 will turn on the display and activate the buttons
 - a. Turning over the dispenser will also do this.
2. OFF: Button 2 will turn off the display and inactivate the buttons

Setup: Carefully Go Through All Steps for Proper Setup

Open the lid: Press your thumb on the front tab, then lift up the cover
CLOCK: Use either a 12-hour clock (AM/PM) or 24-hour clock and set current time. The dispenser's default uses the 24-hour clock

1. Press and hold Button 2 and Button 3 simultaneously to switch to 12-hour mode. "AM" or "PM" will be displayed.
2. Press and hold Button 1 until the hours blink
3. Press Button 2 (reverse) or 3 (forward) to set the correct hours
4. Press Button 1 until the minutes blink
5. Press Button 2 (reverse) or 3 (forward) to set the correct minutes
6. Press and hold Button 1 for 5 seconds to save

ALARM TIME: Alarms are to be set in order, set the earliest first

1. Press Button 1. The hours and "A1" will begin to flash
 - a. Press Button 2 or 3 to set the alarm's hour
 - b. Press Button 1 until the minutes blink
 - c. Press Button 2 or 3 to set the alarm's minute
 - d. Press Button 1 to save and to move to the next alarm
 - e. Repeat a-d until all your alarms are set

2. Delete an Alarm: At the alarm you want to delete, change the hour until --: -- is displayed. Occurs between 11pm & 12am
3. Press 1 until the current time is displayed

ALARM TONE: Tones A, B, C and None (blank) are available. Set Tones with the current time displayed.

1. Press and hold both Button 1 & 3 for a few seconds
2. After the tone changes let go of both buttons completely
3. Repeat 1 & 2 to rotate through the different tones
 - a. After tone C, no text is shown, indicating a light only alarm

SETUP - STEP 2

Filling the Tray:

1. With the lid still open, remove the tray that is inside the dispenser. This makes filling the tray easier.
2. e-pill MedTime Station comes with six different dosing templates. The template identifies the number of alarms per day and, very importantly, the time of the day that you are filling the tray.
3. Fill compartments counterclockwise (up to 28 doses are available) with pills using the dosing template as a guide. Note: If you take medications three (3) times per day, the last compartment will be left empty.

Begin to Use:

1. Align the tray in the dispenser: The next dose will be the first compartment to the right of the compartment placed directly at 6 o'clock.
2. Close the lid
3. Turn the dispenser over and lock the dispenser.
 - a. **WARNING:** When locking/unlocking the cover do not turn the key more than $\frac{1}{4}$ turn.

Congratulations:

All of SETUP – STEP 1 and all of SETUP – STEP 2 is complete. Your e-pill MedTime Station is now 100% setup.

General Information


Stop the Alarm:

Pull the blue handle towards you to release the correct dose and silence the alarm.

Lock:

The dispenser's lock should only be turned $\frac{1}{4}$ of a turn. Overturning the key with a heavy hand can cause the lock to break.

Low Battery Alert:

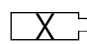 Indicates the batteries are running low. Replace immediately.

Battery Replacement:


The dispenser's clock and alarm times are saved when replacing the batteries, but this is only true if replacement is done prior to the batteries reaching 100% depletion. The new batteries should be inserted within 5 minutes of the old batteries being removed so the dispenser's clock and alarm times stay saved. After replacing the batteries, do check the clock and alarm times to verify their accuracy.

Battery Save:

The dispenser's display will turn off after 5 minutes of inactivity in order to increase battery life. Your medication(s) will still be dispensed at alarm time. The display can be forced to turn on by tipping the dispenser upside down then right side up.

Cleaning:

Remove the medication tray from the dispenser and wipe down the tray with a damp cloth.

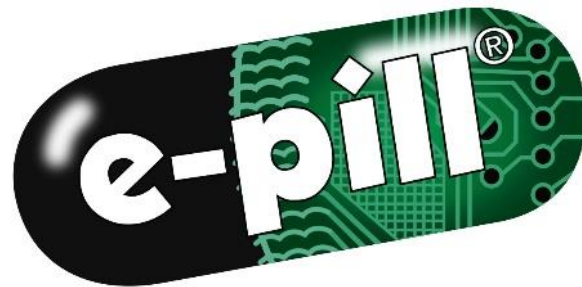

Medication Reminders

Follow us on:



e-pill, LLC
49 Walnut Street, Bldg.4
Wellesley, MA 02481, USA

United States: 1-800-549-0095
International: +1-781-239-2941
www.epill.com