

e-pill[®]

7 Alarm Vibrating Pillbox

SKU 793405



Care and Use Guide

7 Alarm Pill Reminder
Helps simplify complex dosing schedules.




Congratulations on the purchase of your new e-pill® 7 Alarm Vibrating Pillbox (SKU 793405). Please familiarize yourself with its' many useful features and read through the following instructions before you begin to use it. This product is for adult use only. Keep away from children.

1. Viewing lens
2. **SET ALARM**
3. **NUMBER**
4. Compartment

5. **MINUTES**
6. Release Latch
7. **BACK LITE**
8. **HOURS**

9. **VitaCarry®**
10. **DISPLAY**
11. **SET TIME**

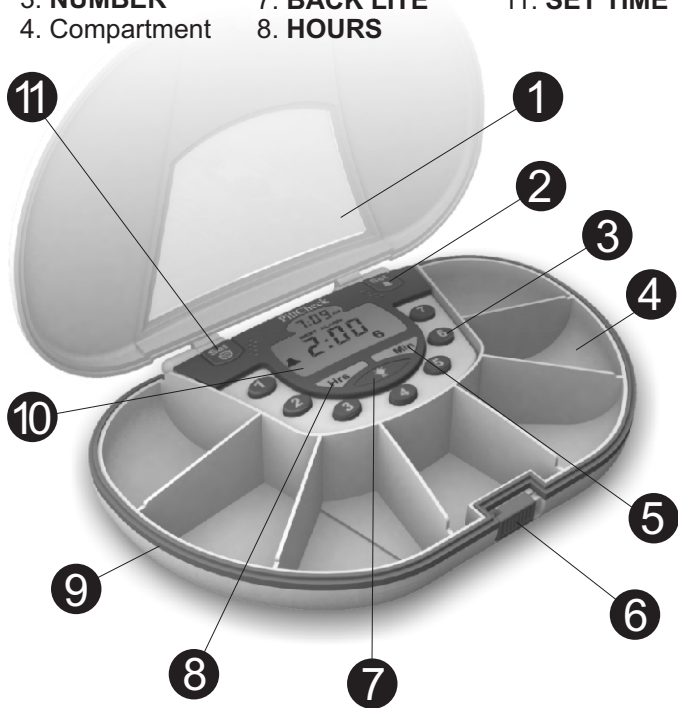

Display at a glance

Time

It is 8:32 AM.



Missed Pill

Compartment 1 was not acknowledged on time and is displayed as a Missed Pill.

Next Alarm

At 10:00 AM the Vibrate and Audible Alarm will activate and remind me to take pills from compartment 3.

Acknowledge

Pressing the number 3 turns off the alarm and the device displays the Next Alarm.(if any)

Before you begin

Read All Instructions








This item is intended for adult use only. Keep medications/supplements away from children. Do not use near water. Do not use for other than the intended purpose. Save these instructions.

Getting started

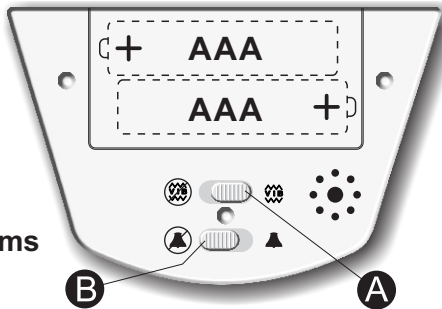Install **AAA Batteries** in back of Timer as shown. You are now ready select initial settings.

Setting Vibrate and Audio Alarms



1. Slide (A) to  or  to turn Vibrate alarm on or off.
2. Slide (B) to  or  to turn audible alarm on or off.

To Set Time

1. Press **(A)** and ^{SET}TIME will appear on the display.
2. Choose **(B)** Hrs & Min.
3. Press **(A)** again to confirm.

Clock

If the AM indicator is *not* displayed, the Time is PM.



Back Lite

Press to enhance display for up to 16 seconds.

To Set Alarm

1. Press **(C)** and ^{SET}ALARM will appear on the display. The alarms 1-7 that are enabled will display and the lowest numbered alarm will begin blinking.
2. For each alarm you wish to set, press a Number Key and set the Alarm Time by using **(B)** Hrs & Min keys.
3. Press **(C)** to confirm your Alarms and exit Set Alarm mode.

To Use Alarm

When a programmed Alarm goes off, it will blink. The Sound and Vibrator may also go on. This is reminding you to take the dose you scheduled for that time and stored in the corresponding compartment. Acknowledge this by pressing the corresponding Number key.

If the Alarm is not acknowledged, Sound and Vibrator will stop after 10 minutes. Alarm will continue to blink. If it is still not acknowledged for up to 1 hour, it will be marked with a X. This indicates you have missed the time you had scheduled to take your dose.

The Next Alarm will be displayed after a Missed Alarm or an acknowledgment. If there are no more Alarms for the day, 0:00 Alarm Time will be displayed.

Helpful hints

Press any Number Key for quick view of that Alarm Time.

Alarm times are in 10 Min intervals to speed up process.

To set Alarm to ZERO, press number twice during Alarm Time set up process.

At midnight, 1st Alarm of the day will be displayed. All Missed pills are reset.

Authorized Reseller

e-pill, LLC

49 Walnut Street, Bldg. 4
Wellesley, MA 02481 USA
1-800-549-0095 | www.epill.com

Manufacturer

eNNOVEA, LLC

2030 Dividend Drive
Columbus, OH 43228
1-513-445-5659 | www.ennovea.com

© eNNOVEA, LLC

VitaCarry® is a US registered trademark.
Designed in USA. Made in China.



Tested to comply
with FCC standards

Manufacturer's Limited Warranty

eNNOVEA LLC will repair or replace any product that fails to operate as promised as a result of defective materials or assembly for a period of one year. Damages caused by mis-use, negligence or accidents are not covered under this warranty. Item should be returned to eNNOVEA, LLC with proof of purchase, postage prepaid and description of problem.