


# OPERATING INSTRUCTIONS

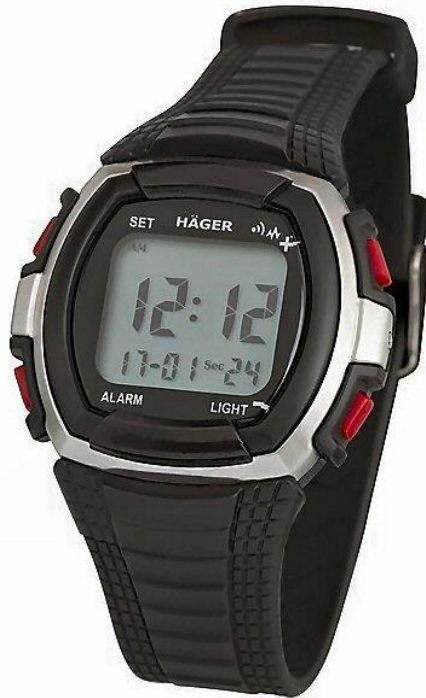

## **e-pill<sup>®</sup> HÄGER<sup>®</sup>**

VIBRATING Alarm / SOUND Alarm / BOTH  
EIGHT (8) ALARM  
MEDICATION REMINDER WATCH  
SKU 931229

Manufactured by HÄGER. Watch Made in China

**e-pill** Medication Reminders

e-pill, LLC | 49 Walnut Street Bldg 4 | Wellesley, MA 02481, USA

Phone (781) 239-2941, Fax (781) 235-3252

Toll-Free 1-800-549-0095

Order your e-pill HÄGER SKU 931229 on-line: [www.epill.com](http://www.epill.com)



# OPERATING INSTRUCTIONS

## e-pill HÄGER Vibrating 8-Alarm Medication Reminder Watch

### BASIC FUNCTIONS

#### FEATURES

The e-pill HÄGER Vibrating Medication Reminder Eight (8) Alarm Watch is an innovative device to help improve medication adherence and patient compliance.

- Up to 8 Daily Alarms
- Choice of Vibration alarm / Sound Alarm or both
- No alarm reset required
- Long alarm duration 3 x 20 seconds (60 seconds in total)
- Easy to read large display
- Back light

#### HOW TO SET THE CORRECT TIME / DATE

Your new e-pill HAGER watch may be in power saving mode (no display). Press A+B+C+D at the same time to start the watch if you do not have a display.

If your watch does not show the correct time please follow these instructions:

- Press & hold button D (SET) for 3 seconds to enter setting mode. The hours will start flashing
- Tap A (+) until the correct Hour is reached. Please note the AM/PM indicator in the top left hand corner. To go backwards tap B (-)
- Tap D
- Tap A until the correct Minute is reached
- Tap D
- Tap A until the correct Year is reached
- Tap D
- Tap A until the correct Month is reached (1=January, 2=February ... 12= December)
- Tap D
- Tap A until the correct Day of the Month is reached
- Tap D to leave setting mode

Congratulations - you have now set your new e-pill HAGER Vibrating 8 Alarm Watch to the correct time and date.

#### TEST MODE

While in normal time display, press & hold button B (LIGHT) for 5 seconds and the vibrator and beeper will turn on for 2 seconds.

#### BACK LIGHT

Tap button B (LIGHT) for Back Light (3 seconds illumination).



#### AM/PM (12 HOUR) or 24 HOUR MODE

Press and hold A button to toggle between AM/PM (12 Hour) or 24 Hour (Military Time) mode.

#### HOW TO SET DAILY ALARMS

While in normal time mode, Tap button C (ALARM) to go to the alarm you want to set. Alarm-1 is AL-1, Alarm-2 is AL-2, Alarm-3 is AL-3 .... and Alarm-8 is AL-8.

For each DAILY ALARM do the following:

- Press & hold button D (SET) for 3 seconds to enter setting mode. The hours will start flashing.
- Tap A until the correct Hour is reached. Please note the AM/PM indicator above in the top left hand corner. To go backwards tap B (-).
- Tap D
- Tap A until the correct Minute is reached.
- Tap D to leave setting mode
- Tap A to select alarm mode for each alarm. Possible settings are Sound Alarm On, Vibrating Alarm On, BOTH (Sound and Vibration) Alarm On and Alarm OFF (no symbol shown)

Proceed to your next alarm - Tap C (ALARM).

Each alarm lasts for 3 x 20 seconds. 20 seconds of alarm (Vibration, Sound or Vibration AND Sound) are followed by 40 seconds of silence. This series repeats three (3) times for a total of 60 seconds. Tap any button to stop the alarm.

#### NO Alarm Reset Required.

The e-pill HAGER will repeat each alarm set on a daily basis – no reset required.


*All the alarms will reset at midnight and automatically repeat every day at the same time.*

## POWER SAVING FUNCTION

You may want to temporarily shutdown your e-pill HÄGER in order to save the battery (power saving mode). Press all four buttons A+B+C+D at the same time while in normal time display. While in power saving mode, the time clock will stop and the display will turn off. Press A+B+C+D to re-start your e-pill HÄGER Alarm Watch.

## BATTERY REPLACEMENT

- Unscrew the watch bottom cover and open it (1).
- Unscrew the battery holder (2). Remove bracket by sliding it to the right.
- Replace the battery with a new 3 Volt Lithium Button Cell (CR2032). The + side facing you (3).
- Replace the battery bracket. Screw back the battery holder.
- Use a metal tool to short circuit the two RESET contacts -see illustration (5).
- Check the watch display (6) to make sure it is working. If not repeat the reset step.
- Screw back the bottom cover. Make sure it is aligned with the words on the bottom cover facing up.
- Perform the TEST MODE function by holding in B button (9) until the vibrator and beeper are turned on.
- Battery life is estimated to be 4 to 6 months but depends in the number of alarms and how long each alarm is sounding.

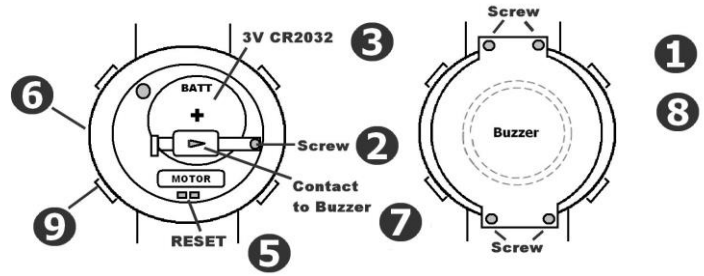

Battery:  
Replaceable CR2032 3V Button Cell

We suggest you take the watch to a jeweler to have the battery replaced.

## WARRANTY

HÄGER covers manufacturing defects for one year from date of purchase. Battery, straps, glass or damage caused by misuse are not covered. In case of defect, send the watch to:

e-pill Medication Reminders  
Legal name: e-pill, LLC  
49 Walnut Street, Building 4  
Wellesley, MA 02481, USA

Or contact your distributor

This watch is designed to help the wearer to remember his/her medication on time. The manufacturer HÄGER takes no responsibility whatsoever for any inconvenience or non-compliance to medications in the event the watch fails to function. The e-pill HÄGER watch is only a supplement to any memory system that you currently use to remind yourself to take your medications or perform medical tasks. The user maintains full responsibility to administer medications or perform medical tasks at the appropriate times. The manufacturer HÄGER or e-pill or the distributor assumes no responsibility for any incidental or consequential damages, including but without limitation, damages or injuries resulting from malfunction of the product, mathematical inaccuracy or loss of alarm times.

### WARNING:

This Operating Instruction is protected by copyright.

©Copyright e-pill, LLC, 2013

e-pill® and HÄGER® are trademarks owned by e-pill, LLC and Monico Properties Ltd.

Exclusive Distributor in the USA, Canada and Mexico

## **e-pill** Medication Reminders

e-pill, LLC, 49 Walnut Street Bldg 4, Wellesley, MA 02481, USA

+1 781 239 2941 or Toll-Free: 1-800-549-0095

Order on-line: [www.epill.com](http://www.epill.com) fax: (781) 235-3252