


# Operating Instructions

## **e-pill** Multi-Alarm Home (759023)

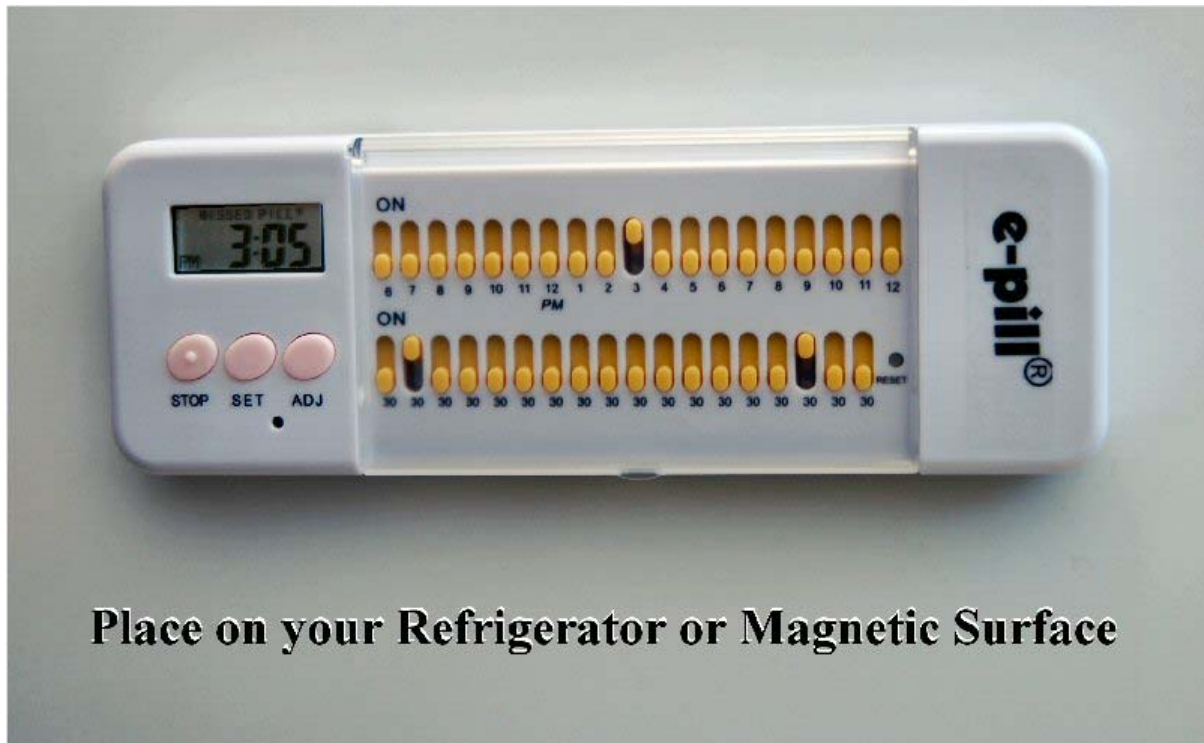

**Place on your Refrigerator or Magnetic Surface**

This e-pill Medication Reminder can remind you to take your medications up to 37 times per day by a beeping alarm. Easy to set alarms repeat automatically every day, no reset required. Alarms can be set from 6 in the morning to midnight on the hour and half hour. Use the magnet to attach this e-pill to your refrigerator door in your kitchen or any magnetic surface in your home.

**Never Forget to take your Medications Again!**

Setting your new e-pill Multi-Alarm HOME:

CLOCK – Set the Current Time

If the clock is not running, please remove the clear insulator strip under the pill tray.

- Press & hold SET button for 3 seconds until the hour start to blink.
- Press the ADJ button to set the hour.
- Press the SET button until the minutes start to blink.
- Press the ADJ button to set minutes.
- Press SET button to finish time setting
- To switch from 12/24-hour clock press the ADJ key after setting the time.

Please turn page

### ALARMS – Set Alarms (up to 37 daily alarms)

In the picture on the front page, the e-pill Multi-Alarm HOME has been set to alarm at the following times: 7:30AM, 3:00PM and 9:30PM. The example to the right shows how to set the alarm for 7:30AM.



# Example: Alarm set for 7:30

- Open the clear cover
- For each alarm, slide ONE switch to the ON (up) position corresponding with your medication time.
- For half hour settings, slide the corresponding half hour(“30”) switch to the ON (up) position. Use the lower level switch only, NOT the top level switch.
- Close the clear cover. Alarms will automatically repeat at the same time each day.

### MEDICATION TIME

#### BEEPING ALARM – What to do when the ALARM goes OFF

- TAKE your MEDICATIONS
- Press STOP to silence alarm after you have taken your medications.

The alarm will sound on and off for a total of 4 minutes. If the STOP button is not pressed within 4 minutes the “MISSED PILL?” indicator will be turned on (flashing on the display).

#### MISSED PILL? Indicator Flashing

If the “MISSED PILL?” indicator is flashing:

- TAKE your MEDICATIONS
- Press STOP to turn off the “MISSED PILL?” indicator.

#### BATTERY – How to replace the Battery

Your e-pill Multi-Alarm Home runs on 2 replaceable AAA batteries. Normal battery life is 12 months. Replace battery when you hear a single beep alarm once every 15 minutes indicating low battery voltage. The low battery voltage indicator icon will also appear on the LCD display. Turn the device over. Unscrew the two screws holding the battery compartment door. Remove battery compartment and old batteries. Insert two fresh AAA alkaline batteries. Put the battery compartment door back on and tighten the screws.

- Press the reset button (gray button under clear cover) and set the correct time.

Disclaimer: This e-pill Medication Reminder is only a supplement to any memory system that you currently use to remind yourself medications. The user maintains full responsibility to administer medication at the appropriate times. e-pill, LLC assumes no responsibility incidental or consequential damages, including but without limitation, damages or injuries resulting from malfunction of the product. e-pill, LLC's liability shall in no circumstances exceed the amounts paid by the user for the timer.



**e-pill**® | Medication Reminders

e-pill, LLC, 70 Walnut Street, Wellesley, MA 02481, USA  
Phone: 1-800-549-0095 (781-239-8255 U.S & International)

[www.epill.com](http://www.epill.com)

Version 060607