



## **e-pill**® Medication Reminders

Distributor:

e-pill, LLC | Medication Reminders  
70 Walnut Street  
Wellesley, MA 02481, USA  
781-239-8255 or Toll-Free: 800-549-0095  
[www.epill.com](http://www.epill.com)  
fax 781-235-3252

Re-order information:

e-pill Vibrating 12 Alarm Pager SKU: 771085  
includes pager, belt clip, leather pouch, manual and batteries



## **OPERATING INSTRUCTIONS**

# **e-pill<sup>®</sup> Vibrating 12 Alarm Pager MEDICATION REMINDER**



## e-pill Vibrating 12 Alarm Pager General Information

Thank you for purchasing an e-pill Medication Reminder. This section provides important general information about your new e-pill Vibrating 12 Alarm Pager.

In this manual we will refer to your e-pill Medication Reminder device as the “Vibrating Pager Reminder”.

### Beep or Vibration Switch (B O V)

Your Vibrating Pager Reminder has a switch on the side labeled **B O V**. This switch is used to make the Vibrating Pager Reminder **Vibrate (V)**, **Beep (B)** or **turn off** the alerts Off (O). The Vibrating Pager Reminder cannot beep and vibrate at the same time (see picture).



### Timer and Alarms - Your Vibrating Pager Reminder has many timing modes

**Press the MODE (red) key to toggle between timing modes**

**12 Alarms (ALARMS SET)** – You can set up to 12 different alarms to vibrate or beep at various times during the day. Use the ENTER button to scroll through set alarms – alarms are displayed in order of next to signal. Any alarm programmed in this mode will signal no matter what other mode the unit is in. No alarm reset is required—the alarms will repeat every day. This is the mode patients find most helpful

**Count Down Timer (TIMER SET)** – (This option is often used by Parkinson’s patients) You can set it to alert once, or set it to repeat itself, vibrating or beeping at any interval you choose for as long as you like. Use for simple timing or for pacing any repetitive need. Great for building new habits, as a memory aid or for pacing presentations.

**Meeting Timer / Interval Timer (MEET-TIMER)** – vibrates at the halfway point, again five minutes before time’s up, and at time’s up, for any time interval you set between 10 – 120 minutes. Each of the three vibrations is distinct. Once time’s up, the *Meeting Timer* automatically resets and repeats on the half hour or hour, depending on the setting. Great for meetings, counseling sessions or presentations.

**Custom Timer (CUSTOM TMR)** – Use to set any countdown time of up to 60 minutes, and program with up to six custom alerts to go off at any time in this period. You can also set it to repeat, or count-up after countdown is complete.

**Stopwatch (STOP-WATCH)** – The *Stopwatch* can be used for counting time for any event. No alarm can be set with the *Stopwatch*, however, any alarm that has been set in the 12 Alarms mode will signal if the unit is in the *Stopwatch* mode.

### Adjustments

**Time and Date adjustment (TIME & DATE)** is used to set the time, date and year.

**Vibration and Beep adjustment (VIB & BEEP)** is used to set the vibration intensity, the beep volume and number of beeps.

### Limited Warranty

The Manufacturer warrants this timer is free from defects in material and workmanship that result in product failure during normal use, according to the following terms and conditions: The limited warranty for this product extends for 1 year beginning on the date of purchase of the product. The limited warranty extends to the original purchaser of the product and is not assignable or transferable to any subsequent purchaser/user. During the limited warranty, TNC will repair, or replace the timer unit, at TNC’s option; or replace any defective parts, or any parts that will not properly operate for their intended use with new items if such repair or replacement is needed because of product malfunction or failure during normal usage. No charge will be made to the consumer for any such parts or repair labor. The limited warranty does not cover defects in appearance, cosmetic, decorative or structural items. TNC shall not be liable for any other losses or damages. These remedies are the consumer’s exclusive remedies for breach of warranty. The consumer shall have no coverage or benefits under this limited warranty if any of the following conditions are applicable: (a) The product has been subject to abnormal use, abnormal conditions, improper storage, exposure to moisture or dampness, unauthorized modifications, unauthorized repair, misuse, neglect, abuse, accident, alteration, improper installation, or other acts which are not the fault of TNC, including damage caused by shipping; (b) the product has been damaged from external causes such as collision with an object, or from fire, flooding, sand, dirt, wind-storm, lightning, earthquake or damage from exposure to weather conditions, an Act of God, or battery leakage, or theft. If a problem develops during the limited warranty period, the consumer shall take the following procedure: Return the product to TNC, 207 Terrace Ave., San Rafael, CA 94901. Include a return address, phone #, description of the problem. TNC will repair or replace the product under the limited warranty within 15 days after receipt of the product by TNC. TNC shall not be liable for special, incidental or consequential damages, including but not limited to, loss of anticipated benefits or profits, loss of savings or revenue, punitive damages, loss of use of the product or any associated equipment, cost of capital, cost of any substitute equipment, downtime, the claims of any third parties including customers, and injury to property, resulting from the purchase or use of the product or arising from breach of the warranty, breach of contract, negligence, strict tort, or any other legal or equitable theory, even if TNC knew of the likelihood of such damages. TNC shall not be liable for delay in rendering service under the limited warranty, or loss of use during the period that the product is being repaired.

Manufacturer:

TNC, 207 Terrace Avenue, San Rafael CA 94901, USA









## 1. Time, Date and Daylight Saving Settings

The e-pill Vibrating Pager Reminder will automatically adjust the time for daylight savings time; please see details below on how to turn this feature on or off). (TIP: Hold Set button down to rapidly advance counting).

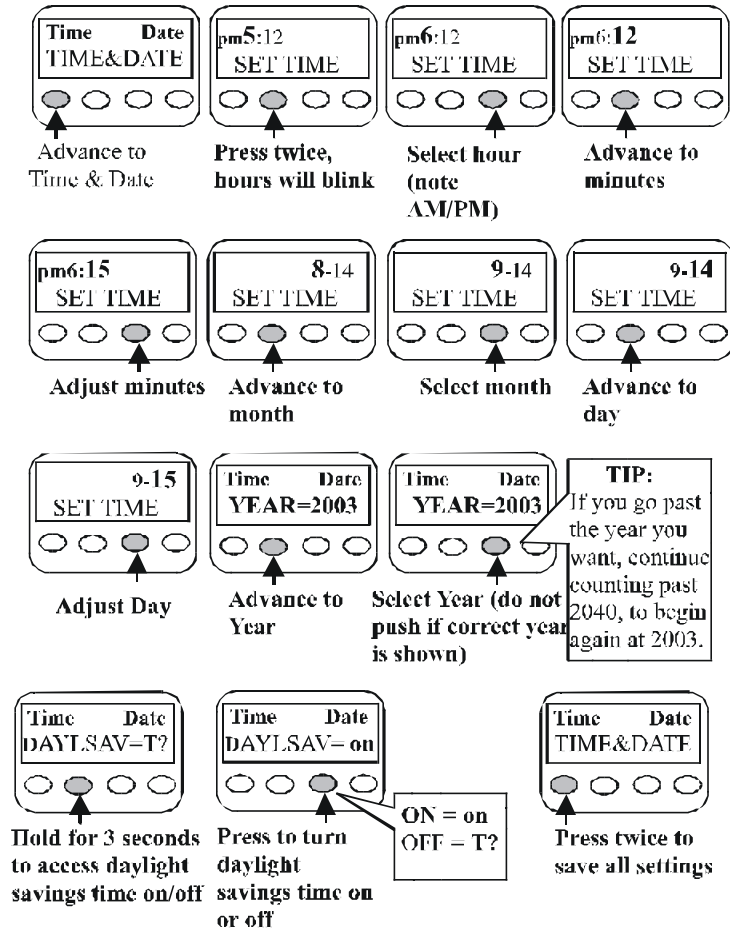

## 2. Vibration Speed, Beep Volume and Number of Beeps

Vibration Speed (1 = low, 2 = medium, 3 = high), Beep Volume (1-3), Number of beeps (1-6). (Factory settings are Vibration = 3; Beep volume = 3; Beep number = 3)

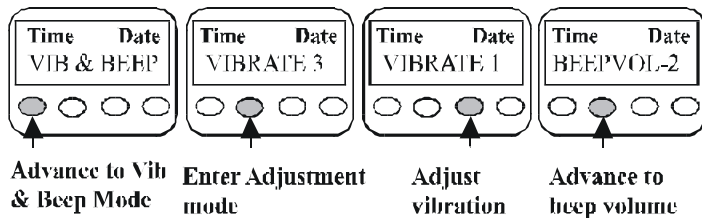

## Q&A

### e-pill Vibrating 12 Alarm Pager

**1. The Unit appears stuck in the Timer or Meeting Timer:** try pushing the two middle buttons simultaneously (just for a second) until the display blinks, and then press the mode button to exit. To EXIT the Custom Timer, press and hold the two middle buttons for 2 seconds.

**2. The Meeting Timer does not appear to be working:** The Meeting Timer does not show countdown when in operation. It only shows "ON" and the total time of the meeting or interval set.

**3. No alert is vibrating for the Meeting Timer:** Make sure you have a meeting time set & the Meeting Timer is ON; the display should show for example: "ON 00:45". Also, make sure the vibration switch on the side of the unit is set to "V" for vibration.

**4. The vibration strength is not strong enough:** Try setting the vibration to "3, high". See "1. How to set vibration speed..." on page 6. You can also try to change battery. Make sure to use an alkaline battery. Note: The e-pill Vibrating Pager will retain your programming for a short time (30 seconds to 1 minute) while you change your battery. See "Care of your Vibrating Pager..." on page 3.

**5. Alarms do not appear to organize correctly when setting multiple alarms:** The unit automatically organizes the set alarms in the order of next to go off. So when you set or check your alarms, the unit will display Alarm 01 as the next to signal, and Alarm 02 as the next after that one and so on.

**6. You cannot enter the Alarms mode to set alarms:** Hold down the Enter button for 3 seconds to enter the Alarms mode (this protects alarms from accidental erasure).

**7. The vibration is too long for the Timer and Alarms:** Vibrations in these modes = 20 seconds. To stop the vibration, push any button. In the Alarms mode, one push of the "Lite" button stops the alert and also engages the snooze feature, and an alert will signal again in 5 minutes.

**8. Can I set alarms with the stop watch:** The Stopwatch cannot be used to set an alarm.



**9. When do I change the Battery:** The battery is low when the backlight becomes dim and the vibration is noticeably slower / or the beep volume is lower. The e-pill Vibrating 12 Alarm Pager will retain your programming for a short time (30 seconds to 1 minute) while you change your battery.

**MEETING TIMER NOTES:**

- Once a meeting time is started, it will repeat every 1/2 hour for durations of 10-30 minutes, and repeat every hour for durations of 31-60 minutes, and repeat every other hour for durations of 61-120 minutes. For example, program a meeting time of 50 minutes and start it on the hour, it will go for 50 minutes, wait 10 minutes and then begin again on the next hour. If you start this 50 minute meeting time at 10 minutes after the hour, it will go for 50 minutes, wait 10 minutes and begin again at 10 minutes after the next hour.
- Meeting times of 30, 60, 120 minutes will repeat back to back automatically. For all other meeting times, there will be a space before the meeting timer repeats. For example, for a meeting time of 25 minutes, it will go for 25 minutes, wait 5 minutes and begin again. For a meeting time of 40 minutes it will go for 40 minutes, wait 20 minutes and begin again. If you want a meeting time to repeat back to back and it is not 30, 60 or 120 minutes, you will need to restart the meeting timer manually at the end of each meeting time by pressing the two middle buttons simultaneously, and then pressing the **Enter** button to re-start the meeting timer at that time.
- When the meeting timer is in operation, it does not show countdown. It only shows "ON" and the total meeting duration.


**7. Stopwatch**

(TIP: You can start the Stopwatch and then leave this mode and it will continue to count. If you do not need the stopwatch to run while using other functions, stop it before leaving this mode.)

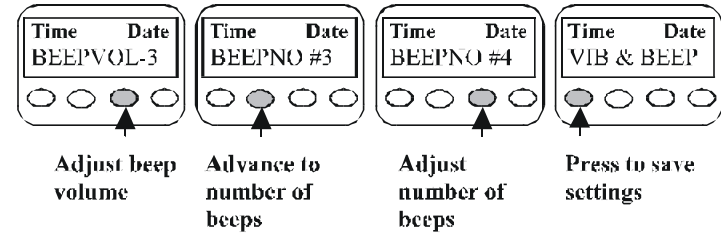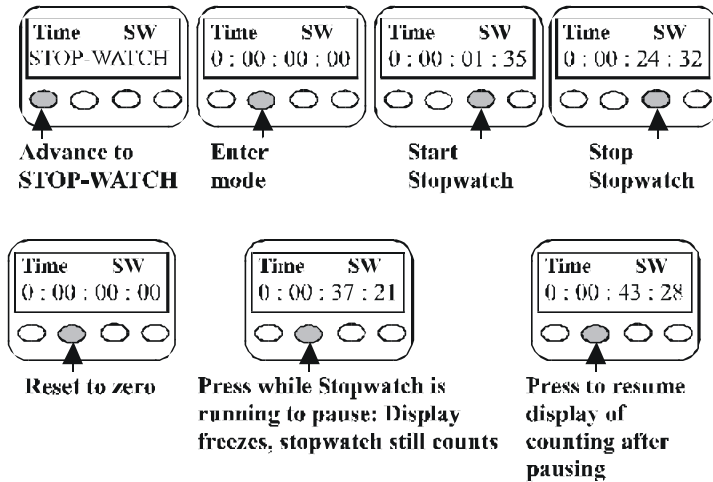

**3. Alarm Settings**

Up to 12 Daily Alarms may be set. **NOTE:** The **Alarm auto-protect** feature requires that you hold the **Enter** button down for 3 seconds to enter the Alarms mode.

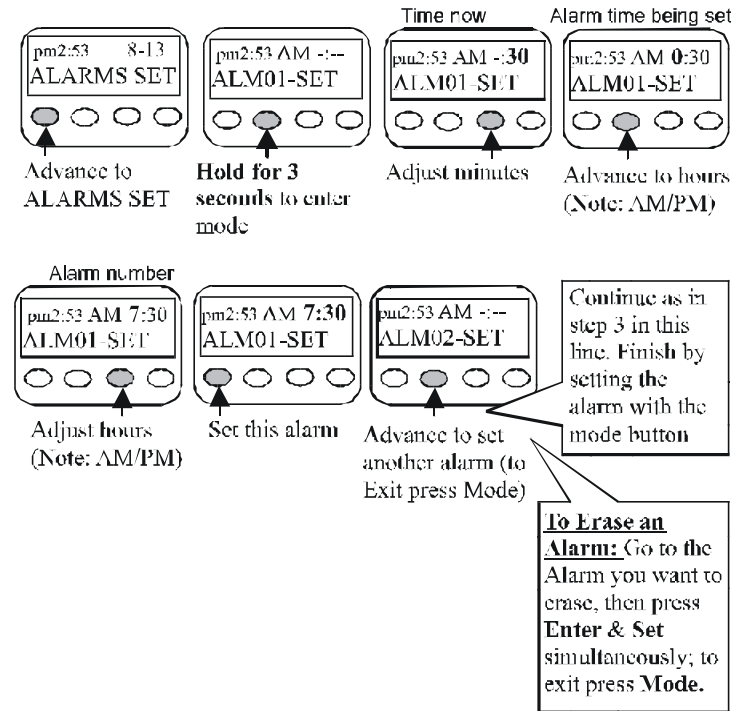

**TIP 1:** To set multiple alarms or to check your alarms, use the **Enter** button to scroll through the alarms. To set additional alarms go to a blank setting, then continue as described above.

**TIP 2:** Vibration alert = 20 seconds. Press any button to stop alert. If you press the **Lite** button, it stops the alert and engages the snooze feature - see **Note 3** below about the snooze button.



**Note 1.** A set alarm will signal in all other modes unless alarm switch (on side of unit) is set to O (off). An alarm remains set after its time has passed, and will sound again every successive day until it is erased.

**Note 2.** Set alarms are automatically organized in the order of next to go off, according to the time of day. In other words, when you check your set alarms, ALARM 01 will be the next to signal, then ALARM 02...etc.


**Note 3. SNOOZE button** - When an alarm signals, press the light button to turn off the alarm and engage the snooze feature; the alarm turns back on automatically after 5 minutes. You can only do this once with each alarm.

**Erasing an Alarm:** To erase an existing alarm, (a) use the **Mode** button to advance to "ALARMS SET"; (b) use the **Enter** button to select the number of the alarm that you wish to erase; (c) Press both the **Enter** and **Set** buttons simultaneously. This will erase the alarm's time, and therefore will erase the alarm itself.


**Changing the Time of an Existing Alarm:** To change the time of an existing alarm, (a) use the **Mode** button to advance to "ALARMS SET"; (b) Use the **Enter** button to select the number of the alarm that you wish to modify; (c) Use the **Set** button to change the minutes, then press the **Enter** button to shift to the hours, which you may change with the **Set** button; (d) Press the **Mode** button to accept the change.

#### 4. Count Down Timer Settings

The Timer can be set to countdown once, or set to repeat over and over.

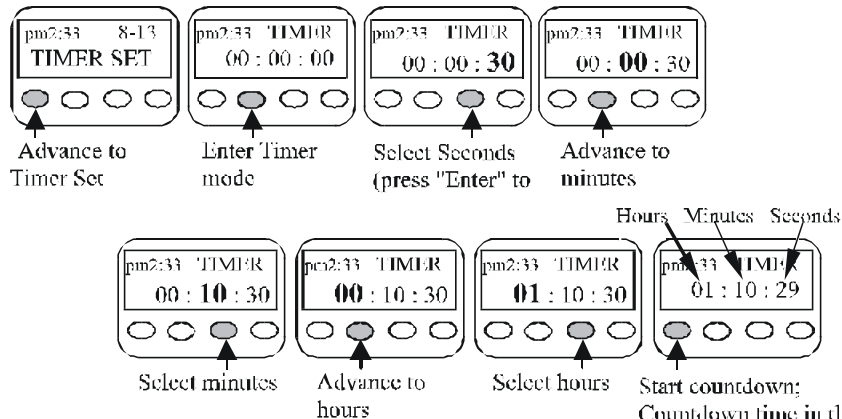


**How to set the Timer to Repeat over & over.**

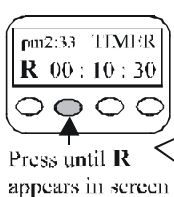

First set the countdown time as shown in the above steps, then before starting the timer, press the "Enter" button until R appears in the left of the display, then press Mode button to start the timer.


**NOTE: Vibration alert lasts for 20 seconds. Press any button to turn off vibration or beep.**

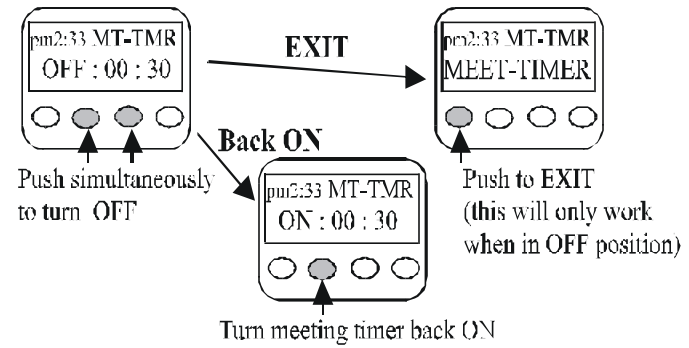

**When the meeting timer is in operation, it does not show countdown. It only shows "ON" and the total meeting duration.**

**Once a meeting time is started, it will repeat every 1/2 hour for durations of 10-30 minutes, and repeat every hour for durations of 31-60 minutes, and repeat every other hour for durations of 61-120 minutes.**


**TIP #1:** The meeting time that was last programmed is saved when you exit. When you return to this mode and wish to re-use this meeting time, follow these steps.

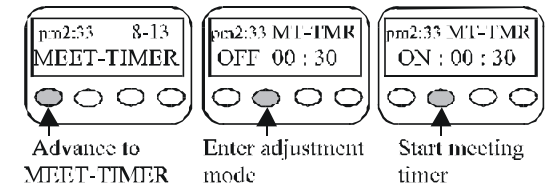

**TIP #2:** To **ERASE** a meeting time go to MEET-TIMER, press **Enter**, press **Set**; this will erase the last used meeting time and allow you to set a new time by pressing **Set** again. Hold **Set** button down for rapid counting.



## Checking Custom Timer programming

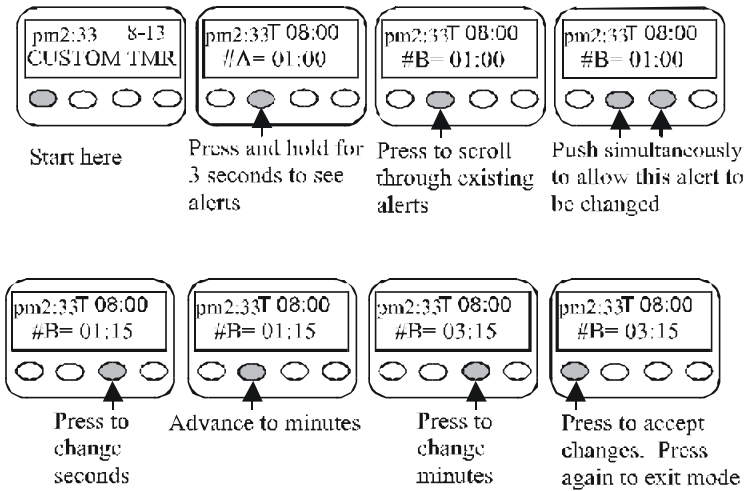


### NOTES for the Custom Timer:


1. When programming, if you have selected all the alerts you need, you can press Mode at any time to save your programming. To start the timer, press **Enter** once, then **Mode** once.
2. Once the custom timer is programmed, it will remain programmed until all programming is erased.
3. To restart timer after countdown completes for one-time countdown or for count-up, press **M** to restart (at any time during count-up).
4. Set any alert to 1 second to have a signal at end of countdown.
5. From "CUSTOM TIMER" you can press **Enter** and wait until ready to start countdown by pressing **Mode**.

## 6. Meeting Timer



Set for any time between 10 minutes and 120 minutes. Signals at 1/2 way through, 5 minutes before the end, and at the end of meeting time. Each signal is unique. Repeats on the 1/2 hour, 1 hour, or every 2 hours.

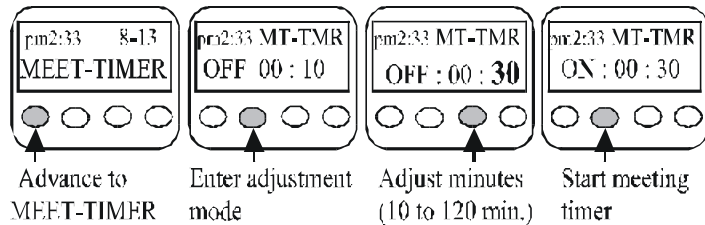

### To EXIT during a countdown

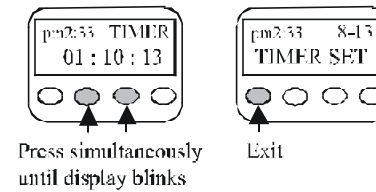


### How to ERASE a countdown time and program a new one

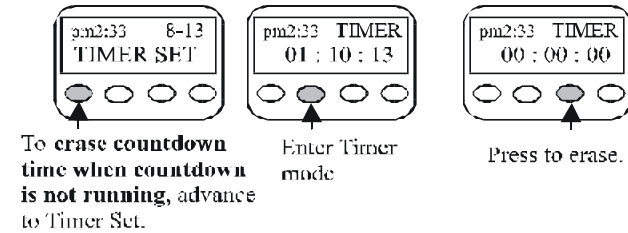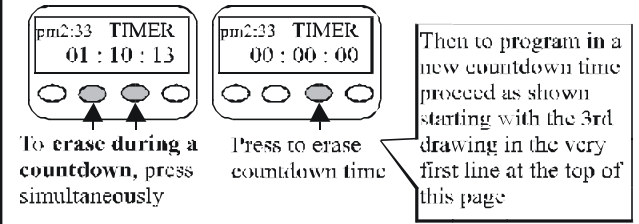


**NOTE:** The Timer saves the last programmed time. When countdown is complete, this saved time is displayed. To use this time again, press **Mode** to start the timer with this time. If you wish to use this time again but have gone to another mode, go to TIMER SET via the **Mode** button, press **Enter** once and press **Mode** to start this countdown time (you can wait as long as you wish to press **Mode** to start the countdown).



## 5. Custom Timer

You can program up to 6 unique alerts to signal anywhere within the countdown time. Set a countdown time of from 2 seconds to 60 minutes. Then within this time, program up to 6 unique alerts to signal at any time within the countdown time period. The countdown time period can be set to countdown once, or set to repeat over and over, or set for a count-up after countdown is finished (so you can see how much time has elapsed since countdown completed).

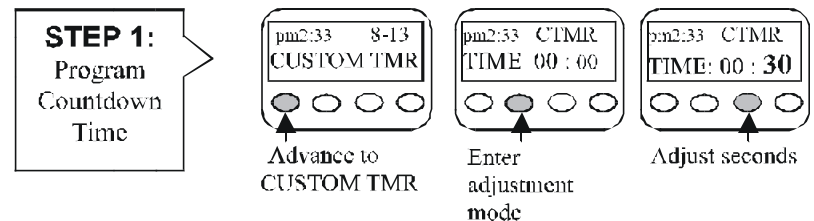

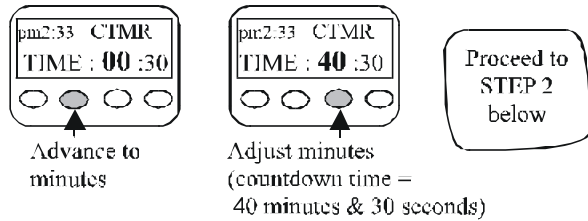


**STEP 2:**  
Program up to 6 alerts within countdown time

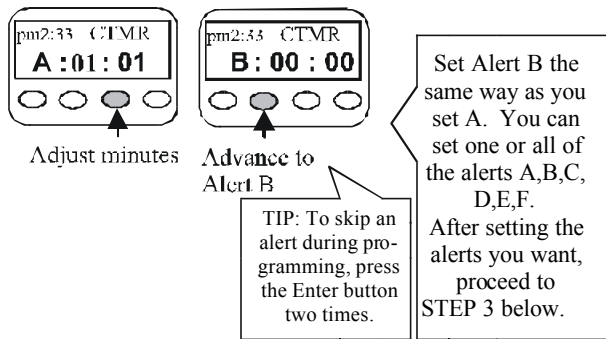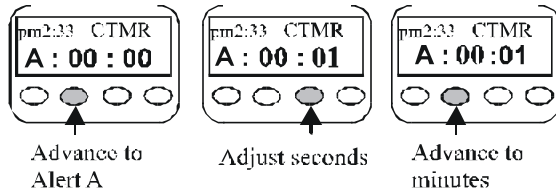


- Alert A = 5 second vibration... or if alert switch (on the side of the unit) is set to beep = 1 beep
- Alert B = 5 second vibration, 5 second space, 5 second vibration.....or..... = 2 beeps
- Alert C = 10 second vibration.....or..... = 3 beeps
- Alert D = 10 second vibration, 7 second space, 10 second vibration.....or..... = 4 beeps
- Alert E = 20 second vibration.....or..... = 5 beeps
- Alert F = 2 second vibration..... or..... = 6 beeps


**STEP 3:** Set to One Countdown, to Repeat, or to Count-Up

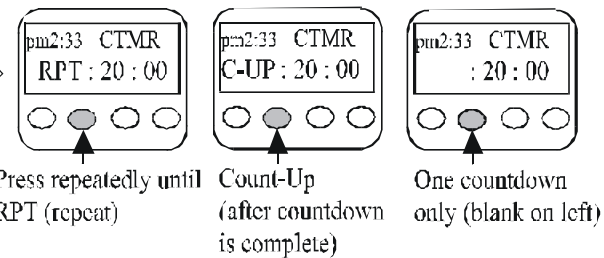


**STEP 4:** Start Countdown (if timer doesn't start, press Enter then press Mode to start)

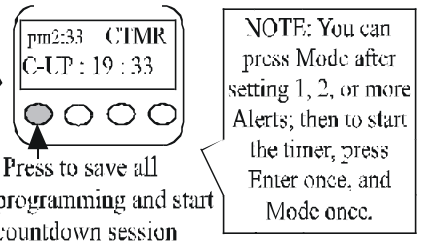


**Stop Countdown**

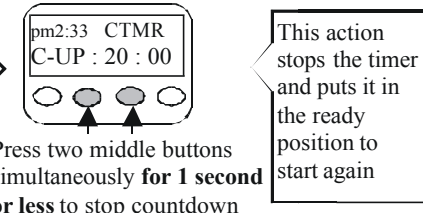


Then **Restart countdown**

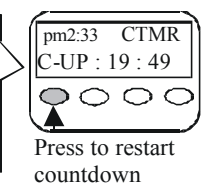

**Exit Mode**

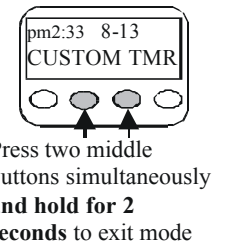

**Restarting Timer when it is programmed**

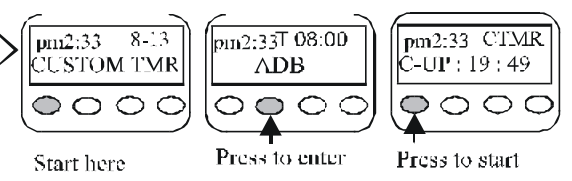

**ERASE ALL programming**

